# *Jota*

Rotierende Dentalinstrumente

Rotary Dental Instruments

Instruments Dentaires Rotatifs

Instrumental Dental Rotatorio

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Josef Rauch, CEO

## Dear Customers and Partners

Over 100 years of experience and the focus on rotary instruments only are a solid based for our offer of high quality products. We align ourselves at your requirements, add innovation and strive after a maximum of customer satisfaction. Our pride and joy is your loyalty, with many of you since decades. Partnership has its responsibility, and we are committed to continue with all our efforts to provide an outstanding market offer of rotary instruments and services. The current catalogue offers new products, extended and completed product lines and a range of useful information around the optimal applications of rotary instruments. You are a valued customer!

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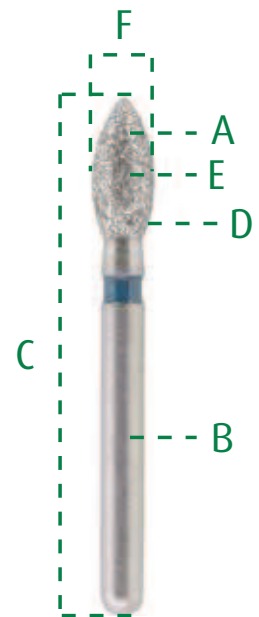
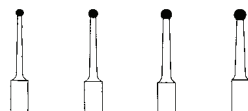
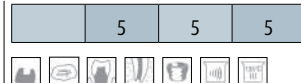
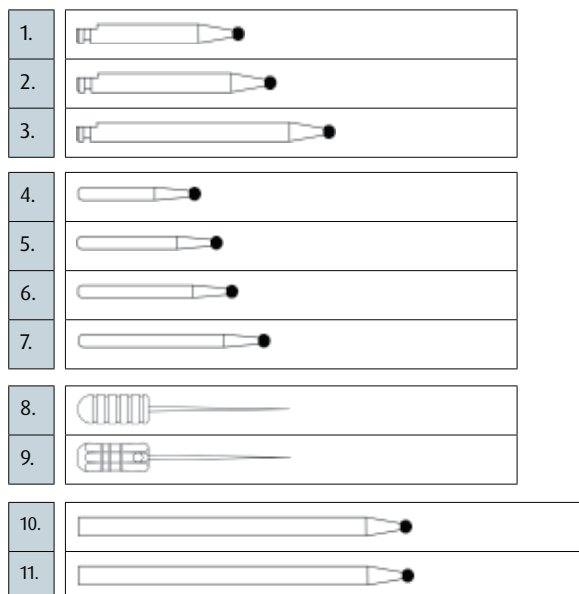
The total lengths of instruments can be longer or shorter according to type of construction.

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› The instruments are protected by practical blister packaging with single tear-off.  
› All information on on figure, size, ISO number and recommended speed are clearly shown.

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The symbols give merely suggestions for the possible implementation of the products. The user decides and takes full responsibility about the precise deployment according to existing indications. Please follow general application and safety instructions for Jota products in the medical and dental area. Details can be found on page 135 as well as on the internet under [www.jota.ch](http://www.jota.ch) or you can request one by mail.

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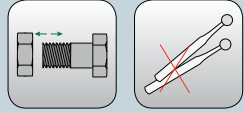
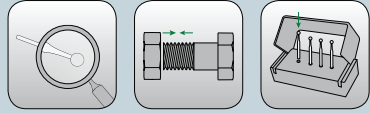
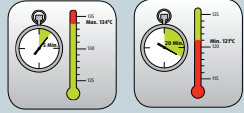
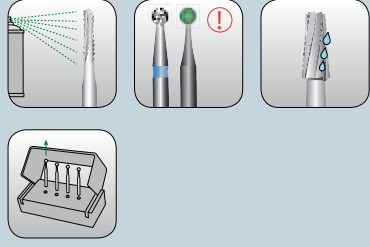
## Manual and Machine cleaning procedure of medical products

Newly delivered unsterile instruments must be prepared accordingly before first use!

1. Pre-treatment	
<ul style="list-style-type: none"> <li>› Abrasive impurities need to be removed from the products directly after use (within two hours maximum). To do so, use running water or a disinfectant solution (must not contain aldehydes, its effectiveness should be established, suitable and compatible for the products)</li> </ul>	
2a. Cleaning and disinfecting	
<ul style="list-style-type: none"> <li>› Disassemble the products</li> <li>› Place the disassembled products in the cleaning bath for the prescribed contact time and concentration, take care that the products do not touch each other</li> <li>› Products with lumen: rinse all instrument lumens 5 times at the beginning and or at the end of the contact time using a disposal syringe and a cannula</li> <li>› Take the products out of the cleaning bath and rinse at least 3 times thoroughly with water</li> <li>› Check the products</li> <li>› Dry the products by blowing them dry using filtered pressurised air (Use only filtered air for drying and water that is either sterile or low in germs and endotoxins)</li> <li>› Wrap the products</li> </ul>	
2b. Cleaning and disinfecting with disinfectant / RDG (MACHINE CLEANING)	
<ul style="list-style-type: none"> <li>› The effectiveness of the disinfectant has to be certified</li> <li>› Disassemble the products as far as possible</li> <li>› Place the disassembled products in the disinfectant and take care that the products do not touch each other</li> <li>› Start the program</li> <li>› Remove the products from the disinfectant after the program has finished</li> <li>› Check and wrap the products straight after removal if possible</li> <li>› Ensure that the disinfectant is regularly maintained and checked. Use only filtered air for drying and water that is either sterile or low in germs and endotoxins</li> </ul>	
3. Checking / Maintenance / Packaging	
<ul style="list-style-type: none"> <li>› Check the products for corrosion, damaged surfaces, bare patches, broken, chipped-off edges, deformations or impurities and eliminate damaged products. Products that are still contaminated need to be cleaned and disinfected once more</li> <li>› Re-assemble disassembled products, instrument oil must not be used</li> <li>› Arrange the cleaned and disinfected products in the dedicated bur block / sterilisation tray.</li> <li>› Wrap the products, bur blocks, sterilisation trays using disposable sterilisation packaging (disposable or double packaging) or sterilisation containers</li> </ul>	
4. Sterilisation	
<ul style="list-style-type: none"> <li>› Steam sterilisation: use of a fractional vacuum process or a gravitation process</li> <li>› Maximum sterilisation temperature 134°C (273°F)</li> <li>› At least 3 min (or 18 min at prion deactivation) at 132°C (270°F) / 134°C (273°F) and drying time of 20 min</li> </ul>	
5. Storage	
<ul style="list-style-type: none"> <li>› The products need to be stored in sterilisation wrapping in a dry and dust-free place. Please note the storage-life resulting from the validation of the sterilisation wrapping</li> </ul>	
i	
<ul style="list-style-type: none"> <li>› New unsteril supplied instruments must be prepared before first use.</li> <li>› Instruments made of tool steel are not suitable for sterilisation without an appropriate pretreatment.</li> <li>› When cleaning instruments unprotected against corrosion (like steel instruments) a corrosion-inhibiting disinfectants and cleaning agents must be used and they have to be pretreated with rust preventing spray before sterilisation.</li> <li>› Furthermore, make sure that instruments of different materials should never be reprocessed together.</li> <li>› Use particular care during the cleaning of the grinding surfaces of the instruments and ensure that all residues are removed by using a brush.</li> <li>› In case of severe contamination of the instrument, it is recommended to use an ultrasonic bath for cleaning.</li> <li>› Usage of protection gloves during work with contaminated instruments is highly recommended!</li> <li>› Bur blocks / instrument trays: <ul style="list-style-type: none"> <li>› Cleaning and disinfecting only without products being loaded (products must not be cleaned and disinfected whilst they are in the bur block / instrument tray)</li> <li>› Bur blocks and other instruments made of aluminum (or other light metal) are not suitable for cleaning using a disinfectant RDG</li> </ul> </li> <li>› Please follow the manufacturer's instructions and the maintenance specifications for using the autoclave and the disinfectant.</li> <li>› The method of use, reaction time and suitability of disinfectants and cleaning agents for certain types of instrument are covered by the manufacturers' instructions.</li> </ul>	
<p>These is a comprehensive discription of our detailed reprocessing advices. This can be found on our website, in our catalogue on page 135 and on request we'd be pleased to send you the printed information.</p>	

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The recommended speeds are a theoretical value, you can find the specific recommended speeds in the product-related information.

To over exceed the given rotation speed can lead to impairment of the efficiency of the products and/or be of danger for the patient as well as the user. The choice of the working speed for a specific application case is determined by the material subsequently to be used, the instrument power, and the maximal down force pressure ultimately decided upon by the user.

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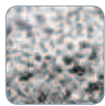
The use of coarse-grained diamonds (ISO 534, 544 and 554) can lead to increased heat generation. When using these products, special care should therefore be taken to ensure sufficient cooling and minimal application force. Instruments from ISO size 031 should be used with additional cooling (water spray)

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## Types and shapes of diamond instruments for dentistry

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- › one-piece tungsten carbide instruments
- › special cutting edge geometry
- › Greatly reduced preparation and treatment times
- › CBP instruments correspond to the most popular diamond-coated preparation instruments for shoulder, chamfer and lingual preparations
- › RRC31R ideal for use as a crown cutter for metal ceramic crowns

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Types and shapes of steel instruments for dental surgeries

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These instruments should be conditioned using a rust protection spray to avoid the formation of black spots/margins or oxidation areas during preparation of non-stainless steel instruments in an autoclave/ steam steriliser.

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### Types and shapes of Abrasives for dental surgeries

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Contact pressure:  
 › recom.: approx. 2N(approx. 200g)

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Contre-angles RA:  
 › recom.: 15'000 - 25'000 min<sup>-1</sup>  
 › max.: 30'000 min<sup>-1</sup>

Turbine FG:  
 › recom.: 80'000-100'000 min<sup>-1</sup>  
 › max.: 120'000 min<sup>-1</sup>

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- › For chairside therapy with glass ceramic like IPS e.max CAD
- › For optimally preparing the restoration for glazing. Smooths finishing striations on the surface and simultaneously roughens the surface, so that the glaze does not run during application.
- › Also suitable for trimming bar attachments!
- › Can also be used for adjusting occlusal contact points with fully anatomical zirconia crowns.
- › Used in a blue contra-angle

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silicon carbide for preparing porcelain and acrylics

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high-grade corundum for preparing medium-hard metal alloys and metal

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Finest grinding/polishing of filling materials and all composite-filling-materials

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Types and shapes of polishers and brushes for dental surgeries

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JOTA EASY - One-step high-performance diamond polisher for finishing and mirror finishing of all composite  
Skimming: 7'000 - 10'000 rpm // Polishing: 3'000 - 8'000 rpm

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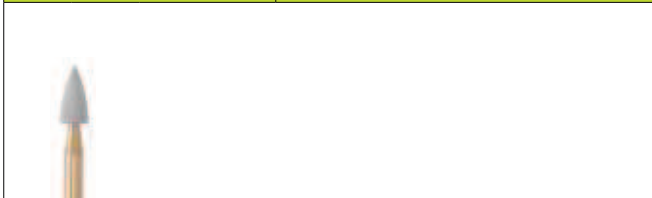
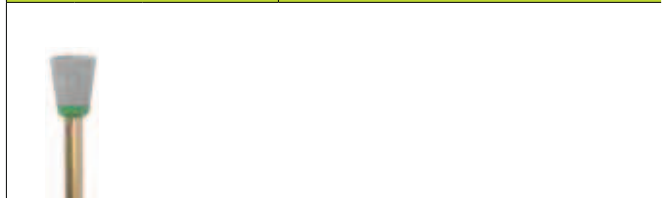
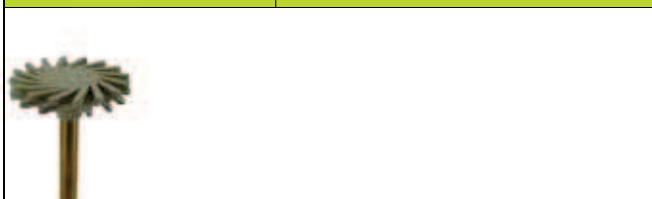
Diamond polisher, which can be optionally used as a preliminary polishing stage prior to using JOTA EASY to produce even more professional results.  
 JOTA PROFESSIONAL Base  
 Step 1 - Green: contouring and smoothing 3'000 - 8'000 rpm  
 JOTA EASY  
 Step 2 - Grey: high-lustre polishing  
 Also available in the kit 1394, see page 129

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**Silicone polishers**

- › two-step silicone polishing system for material smoothing and polishing of composite
- › light blue: smoothing and pre-polishing / yellow: polishing

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**Diamond polishers**

- › three-step diamond polishing system for material removal, smoothing and high-lustre polishing of porcelain
- › also suitable for metals - after use on metals, clean the polishers to prevent striation marks forming on the porcelain
- › blue: material removal / pink: smoothing and pre-polishing / grey: high-lustre polishing

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**Silicone polishers**

- › three-step silicone polishing system for material removal, smoothing and polishing of porcelain
- › light grey: material removal / pink: smoothing and pre-polishing / grey: polishing

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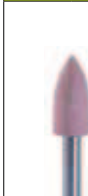
**Diamond polishers**  
 › two-step diamond polishing system for material removal, smoothing and gloss polishing of lithium disilicate  
 › red-brown: smoothing and pre-polishing / beige: gloss polishing

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## Recommendation

Rotate the trephine a few times anti-clockwise to scribe a line. Then cut to the full depth intermittently in a clockwise direction.

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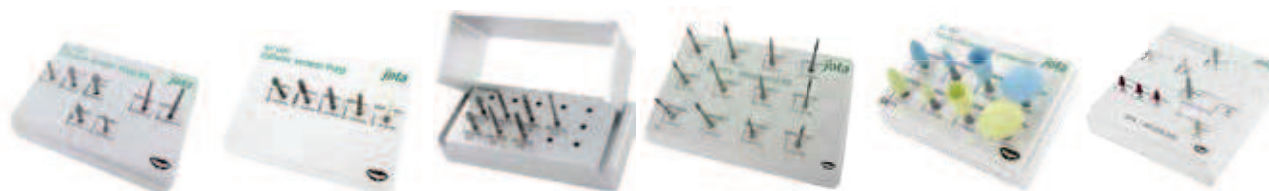
1418

JOTA Kit Crown prep conventional  
JOTA Kit Crown prep conventional

Fig	Shank	ISO	Ø							
<b>Turbine   Friction Grip</b>										
882	FG	806 314 142 524 -	016							
882F	FG	806 314 142 514 -		014						
833	FG	806 314 277 524 -			023					
508G	FG	806 314 508 534 -				020				
508F	FG	806 314 508 514 -				020				
878F	FG	806 314 298 514 -						016		
859	FG	806 314 166 524 -								014
			1	1	1	1	1	1	1	1

# 1424

Universal Prep Kit - University of Bern

Fig	Shank	ISO	Ø			
<b>Turbine   Friction Grip</b>						
833	FGL	806 315 277 524 -				023
833F	FGL	806 315 277 514 -				023
558	FG	806 314 558 524 -		013	015	
558F	FG	806 314 558 514 -		013	015	
859L	FG	806 314 167 524 -	012			
893F	FG	806 314 507 514 -				023
			1	1	1	1

# 1430

 JOTA Kit Inlay/Onlay Adhesive  
 JOTA Kit Inlay/Onlay Adhesive

Fig	Shank	ISO	Ø			
<b>Turbine   Friction Grip</b>						
839	FG	806 314 150 524 -	016			
845RF	FG	806 314 544 514 -		025		
846	FG	806 314 171 524 -			016	
846	FG	806 314 171 524 -			025	
846RF	FG	806 314 545 512 -				016
846REF	FG	806 314 545 504 -				016
852F	FG	806 314 199 514 -				014
			1	1	1	1

# 1431

 JOTA Universal All Ceramic Kit, nach Prof. Brägger, Universität Bern  
 JOTA Universal All Ceramic Kit, acc. Prof. Brägger, University of Bern

Fig	Shank	ISO	Ø			
<b>Turbine   Friction Grip</b>						
852	FG	806 314 199 524 -	014			
852	FG	806 314 199 524 -		016		
859L	FG	806 314 167 524 -			012	
833	FG	806 314 277 524 -			023	
525	FG	806 314 525 524 -				017
526	FG	806 314 526 524 -				023
852F	FG	806 314 199 514 -	014			
852F	FG	806 314 199 514 -		016		
833F	FG	806 314 277 514 -			023	
893F	FG	806 314 507 514 -				023
			1	1	1	1

## 1435

Jota LS Gloss - Dentistry

Fig	Shank	ISO	Ø			
<b>Turbine   Friction Grip</b>						
LS9871M	RA	803 204 243 522 -	040			
LS9872M	RA	803 204 030 522 -		060		
LS9871F	RA	803 204 243 512 -			040	
LS9872F	RA	803 204 030 512 -				060
			1	1	1	1

## 1436

Jota ZIR Gloss - Dentistry

Fig	Shank	ISO	Ø			
<b>Turbine   Friction Grip</b>						
ZIR9861M	RA	803 204 243 521 -	040			
ZIR9862M	RA	803 204 030 521 -		060		
ZIR9861F	RA	803 204 243 511 -			040	
ZIR9862F	RA	803 204 030 511 -				060
			1	1	1	1

## 1439

Nano Ceramic Chairside Polishing

Fig	Shank	ISO	Ø			
<b>Turbine   Friction Grip</b>						
9837	RA	803 204 543 503 -	140			
9150	RA	652 204 546 503 -		140		
GP662	RA	655 204 288 526 -			035	
830LF	RAL	806 205 258 514 -				014
			1	1	1	1

# 1440

## LS GLOSS CHAIRSIDE KIT

Fig	Shank	ISO	Ø					
<b>Turbine   Friction Grip</b>								
LS9873M	RA	803 204 243 522 -	040					
LS9873F	RA	803 204 243 512 -		040				
833F	FG	806 314 277 514 -			023			
GP662	RA	655 204 288 526 -				035		
850F	FG	806 314 198 514 -						016
			1	1	1	1	1	1

# 1441

## ZIR GLOSS CHAIRSIDE SET

Fig	Shank	ISO	Ø					
<b>Turbine   Friction Grip</b>								
ZIR9863M	RA	803 204 243 521 -	040					
ZIR9863F	RA	803 204 243 511 -		040				
GP662	RA	655 204 288 526 -			035			
Z833F	FG	806 314 277 314 -				023		
Z850F	FG	806 314 198 314 -						018
			1	1	1	1	1	1

# 1443

## Efficient Veneer Prep Kit

Fig	Shank	ISO	Ø					
<b>Turbine   Friction Grip</b>								
801	FG	806 314 001 524 -	014					
815	FG	806 314 042 524 -		022				
815	FG	806 314 042 524 -			026			
833F	FG	806 314 277 514 -				023		
890EF	FG	806 314 160 504 -					010	
852GF	FG	856 314 199 534 -						014
852GF	FG	856 314 199 534 -						016
			1	1	1	1	1	1

## 1449

## Esthetic Veneer Prep Kit

Fig	Shank	ISO	Ø			
-----	-------	-----	---	--	--	--

## Turbine | Friction Grip

801	FG	806 314 001 524 -	018				
834	FG	806 314 552 524 -		018			
852	FG	806 314 199 524 -			016		
852F	FG	806 314 199 514 -				016	
649	FG	635 314 171 505 -					025
			1	1	1	1	1

## 1460

JOTA ZirPrep Kit  
JOTA ZirPrep Kit

Fig	Shank	ISO	Ø			
-----	-------	-----	---	--	--	--

## Turbine | Friction Grip

Z801L	FG	806 314 697 324 -	014				
Z833	FG	806 314 277 324 -		023			
Z833F	FG	806 314 277 314 -		023			
Z850	FG	806 314 198 324 -			018		
Z850F	FG	806 314 198 314 -			018		
Z863	FG	806 314 250 324 -					012
Z863F	FG	806 314 250 314 -					012
			1	1	1	1	1

# 1771

## Jota Paradontal Kit

Fig	Shank	ISO	Ø									
<b>Turbine   Friction Grip</b>												
830LF	RAL	806 205 258 514 -	014									
830LEF	RAL	806 205 258 504 -		014								
830LUF	RAL	806 205 258 494 -			014							
830LF	RAXL	806 206 258 514 -				014						
830LEF	RAXL	806 206 258 504 -					014					
830LUF	RAXL	806 206 258 494 -						014				
849LF	RAL	806 205 693 514 -							012			
849LEF	RAL	806 205 693 504 -								012		
849LF	RAXL	806 206 693 514 -									012	
849LEF	RAXL	806 206 693 504 -										012
849LUF	RAXL	806 206 693 494 -										012
849LUF	RAL	806 205 693 494 -										012
			1	1	1	1	1	1				

# 1899

## Jota Kit Silicon Polisher Composite

Fig	Shank	ISO	Ø									
<b>Turbine   Friction Grip</b>												
9121M	RA	685 204 243 522 -	050									
9122M	RA	685 204 243 522 -		030								
9123M	RA	685 204 030 522 -			070							
9124M	RA	685 204 304 522 -				110						
9121F	RA	685 204 243 512 -					050					
9122F	RA	685 204 243 512 -						030				
9123F	RA	685 204 030 512 -							070			
9124F	RA	685 204 304 512 -										110
			1	1	1	1	1	1	1			

# 1912

## Periimplantitis/Parodontitis Workshop

Fig	Shank	ISO	Ø									
<b>Turbine   Friction Grip</b>												
811LG	FG	806 314 039 534 -	037									
C1S	RA	500 204 001 003 -		023								
C1S	RAL	500 205 001 003 -			023							
652	RA	655 204 173 523 -				035						
638	RA	635 204 110 505 -					025					
9501M	RA	653 204 243 513 -							030			
9831	RA	803 204 243 502 -										030
			1	1	1	1	1	1	1			

Zubehör

Accessories

Accessoires

Accesorios

**Quick-n-Clean**
**Quick-n-Clean**
**Quick-n-Clean**
**Quick-n-Clean**

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**The JOTA Quick-n-Clean is a concentrate used for cleaning and disinfection of dental and medical instruments.**
**Properties:**

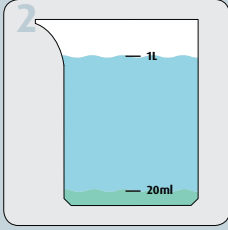
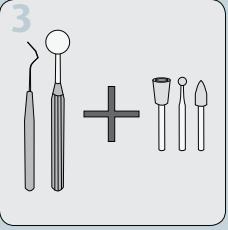
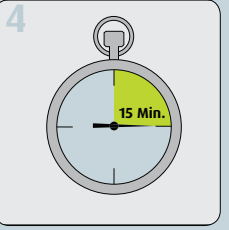
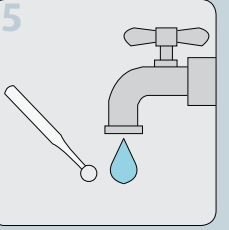
- › 2% concentration
- › Proven effect within 15 minutes against bacteria, fungi and viruses
- › Cost-effective: 1 L concentrate provides 50 L solution
- › Very good material compatibility, even with sensitive materials like steel, aluminum, rubber and silicone
- › Contains highly effective corrosion inhibitor
- › Depending on the application the ready-to-use solution can be used for up to one week, even when diluted
- › Suitable for immersion solution or ultrasonic cleaner
- › DGHM/VAH registered
- › Dermatologically tested
- › Aldehyde-free and phenol-free

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## Instructions for processing &amp; General application and safety instructions

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### Instructions for the processing (cleaning, disinfection, and sterilization) of instruments from Jota AG

Issued: January 2017

The medical devices produced and sold by Jota AG are re-usable unless their label contains explicit information to the contrary. However, as a rule, it is the sole responsibility of the doctor/expert using the devices to decide whether, depending on the respective case and the potential wear and tear of the products, he can re-use the products and how frequently he uses them. In case of doubt, it is always advisable to discard the products early and to replace them. The manufacturer Jota AG cannot guarantee the faultless function and performance of the products combined with a maximum degree of safety if the products are overused. These reprocessing instructions apply in principle to all medical devices making up the product range supplied by Jota AG. Any particular features and/or exclusions that only concern individual items or groups of items are referred to separately.

#### Fundamental points

All instruments are to be cleaned, disinfected, and sterilized prior to each application; this is required as well for the first use after delivery of the unsterile instruments (cleaning and disinfection after removal of the protective packaging, sterilization after packaging). An effective cleaning and disinfection is an indispensable requirement for an effective sterilization of the instruments.

You are responsible for the sterility of the instruments. Therefore, please ensure that only sufficiently device and product specifically validated procedures will be used for cleaning, disinfection, and sterilization, that the used devices (WD, sterilizer) will be maintained and checked regularly, as well as that the validated parameters will be applied for each cycle.

Please pay attention to avoid a higher contamination of the complete bur block during application; otherwise it is necessary to clean and disinfect the bur block as well as all instruments inside (after removal).

Additionally, please pay attention to the legal provisions valid for your country as well as to the hygienic instructions of the doctor's practice or of the hospital. This applies particularly to the different guidelines regarding the inactivation of prions (not relevant for USA).

**Some instruments require additional aspects. For this, pay attention to chapter "Specific aspects".**

#### Cleaning and disinfecting

##### Basic:

If possible, an automated procedure (WD (Washer-Disinfector)) should be used for cleaning and disinfection of the instruments. A manual procedure – even in case of application of an ultrasonic bath – should only be used if an automated procedure is not available; in this case, the significantly lower efficiency and reproducibility of a manual procedure has to be considered.

The pre-treatment step is to be performed in both cases.

##### Pre-treatment:

Please remove coarse impurities of the instruments directly after application (within a maximum of 2 h).

##### Procedure:

1. Rinse the instruments at least 1 min under running water (temperature < 35 °C/95 °F).
2. Soak the instruments at least for the given soaking time in the pre-cleaning solution<sup>1</sup> (by the use of an ultrasonic bath) so that the instruments are sufficiently covered. Pay attention that there is no contact between the instruments. Assist cleaning by careful brushing with a soft brush (at least three times after at beginning of soaking, aids see chapter "Specific aspects").
3. Activate ultrasonic treatment for an additional soaking time (but not less than 5 min).
4. Then, remove the instruments of the pre-cleaning solution and post-rinse them at least three times intensively (at least 1 min) with water.
5. In case of still visible contamination repeat steps 2, 3, and 4, otherwise discard the instrument. This is especially relevant for diamond instruments.

Pay attention to following points during selection of the cleaning detergent:

- › fundamental suitability for the cleaning of instruments made of metallic or plastic material
- › suitability of the cleaning detergent for ultrasonic cleaning (no foam development)
- › compatibility of the cleaning detergent with the instruments (see chapter „material resistance,,)

Pay attention to the instructions of the detergent manufacturer regarding concentration, temperature and soaking time as well as post-rinsing. Please use only freshly prepared solutions as well as only demineralized sterile or low contaminated water (max. 10 germs/ml) as well as low endotoxin contaminated water (max. 0.25 endotoxin units/ml), for example purified/highly purified water, and a soft, clean, and lint-free cloth and/or filtered air for drying, respectively.

<sup>1</sup> In case of application of a cleaning and disinfection detergent for this (e.g. in consequence of personnel's safety) please consider, that this should be aldehyde-free (otherwise fixation of blood impurities), possess a fundamentally approved efficiency (for example VAH/DGHH or FDA/EPA approval/clearance/registration or CE marking), be suitable for the disinfection of instruments made of metallic or plastic material, and be compatible with the instruments (see chapter „material resistance,,). Please consider, that a disinfectant used in the pre-treatment step serves only the personnel's safety, but cannot replace the disinfection step later to be performed after cleaning.

#### Automated cleaning/disinfection (WD (Washer-Disinfector)):

Pay attention to following points during selection of the WD:

- › fundamentally approved efficiency of the WD (for example CE marking according to EN ISO 15883 or DGHH or FDA approval/clearance/registration)
- › possibility for an approved program for thermal disinfection (A0 value > 3000 or – in case of older devices – at least 5 min at 90°C/194 °F; in case of chemical disinfection danger of remnants of the disinfectant on the instruments)
- › fundamental suitability of the program for instruments as well as sufficient rinsing steps in the program
- › post-rinsing only with demineralized sterile or low contaminated water (max. 10 germs/ml, max. 0.25 endotoxin units/ml), for example purified/highly purified water
- › only use of filtered air (oil-free, low contamination with microorganisms and particles) for drying
- › regularly maintenance and check/calibration of the WD

Pay attention to following points during selection of the cleaning detergent:

- › fundamental suitability for the cleaning of instruments made of metallic or plastic material
- › additional application – in case of non-application of a thermal disinfection – of a suitable disinfectant with approved efficiency (for example VAH/DGHH or FDA/EPA approval/clearance/registration or CE marking) compatible to the used cleaning detergent
- › compatibility of the used detergents with the instruments (see chapter „material resistance,,)

Pay attention to the instructions of the detergent manufacturers regarding concentration, temperature and soaking time as well as post-rinsing.

##### Procedure:

1. Transfer the instruments in the WD by the use of a small pieces basket.
2. Start the program.
3. Remove the instruments of the WD after end of the program.
4. Check and pack the instruments immediately after the removal (see chapters „check,, „maintenance,, , and „packaging,, , if necessary after additional post-drying at a clean place).

The fundamental suitability of the instruments for an effective automated cleaning and disinfection was demonstrated by an independent, governmentally accredited and recognized (§ 15 (5) MPG) test laboratory by application of the WD G 7836 CD, Miele & Cie. GmbH & Co., Gütersloh, (thermal disinfection) and the pre-cleaning and cleaning detergent Neodisher mediclean forte (Dr. Weigert GmbH & Co. KG, Hamburg) considering to the specified procedure.

#### Manual cleaning and disinfection:

Pay attention to following points during selection of the cleaning and disinfection detergents:

- › fundamental suitability for the cleaning and disinfection of instruments made of metallic or plastic material
- › in case of application of an ultrasonic bath: suitability of the cleaning detergent for ultrasonic cleaning (no foam development)
- › application of a disinfectant with approved efficiency (for example VAH/DGHH or FDA/EPA approval/clearance/registration or CE marking) compatible with the used cleaning detergent
- › compatibility of the used detergents with the instruments (see chapter „material resistance,,)

Combined cleaning/disinfection detergents should not be used.

Only in case of extremely low contamination (no visible impurities) combined cleaning/disinfection could be used.

Pay attention to the instructions of the detergent manufacturers regarding concentration, temperature and soaking time as well as post-rinsing. Please use only freshly prepared solutions as well as only demineralized sterile or low contaminated water (max. 10 germs/ml) as well as low endotoxin contaminated water (max. 0.25 endotoxin units/ml), for example purified/highly purified water, and a soft, clean, and lint-free cloth and/or filtered air for drying, respectively.

##### Procedure:

###### Cleaning

1. Soak the instruments for the given soaking time in the cleaning solution (by the use of an ultrasonic bath) so that the instruments are sufficiently covered. Pay attention that there is no contact between the instruments. Assist cleaning by careful brushing with a soft brush (at least three times after at beginning of soaking, aids see chapter "Specific aspects").
2. Activate ultrasonic treatment for an additional soaking time (but not less than 15 min).
3. Then, remove the instruments of the cleaning solution and post-rinse them at least three times intensively (at least 1 min) with water.
4. Check the instruments (see chapters „check,, and „maintenance,,).

###### Disinfection

5. Soak the instruments for the given soaking time in the disinfectant solution so that the instruments are sufficiently covered. Pay attention that there is no contact between the instruments.
6. Then, remove the instruments of the disinfectant solution and post-rinse them at least five times intensively (at least 1 min) with water.
7. Dry and pack the instruments immediately after the removal (see chapter „packaging,, , if necessary after additional post-drying at a clean place).

The fundamental suitability of the instruments for an effective cleaning and disinfection was demonstrated by an independent, governmentally accredited and recognized (§ 15 (5) MPG) test laboratory by application of the pre-cleaning and cleaning detergent Cidezyme/Enzol and the disinfectant Cidex OPA (Johnson & Johnson

GmbH, Norderstedt) considering the specified procedure.

##### Check

Check all instruments after cleaning or cleaning/disinfection, respectively, on corrosion, damaged surfaces, and impurities. Do not further use damaged instruments (for limitation of the numbers of re-use cycles see chapter „reusability,,). Still dirty instruments are to be cleaned and disinfected again.

##### Maintenance

Instrument oils or grease must not be use with the exception of steel instruments. In that case use only instrument oils (white oil) admitted to steam sterilization considering the maximum possible sterilization temperature, with approved biocompatibility and without mono-, di, or triethanolamine as corrosion inhibitor.

##### Packaging

Please insert the cleaned and disinfected instruments in the corresponding bur blocks (if required) and pack them in single-use sterilization packagings (single or double packaging), which fulfill the following requirements (material/process):

- › EN ISO/ANSI AAMI ISO 11607 (for USA: FDA clearance)
- › suitable for steam sterilization (temperature resistance up to at least 142 °C (288 °F), sufficient steam permeability)
- › sufficient protection of the instruments as well as of the sterilization packagings to mechanical damage

##### Sterilization

Please use for sterilization only the listed sterilization procedures; other sterilization procedures must not be applied.

##### Steam sterilization

- › fractionated vacuum/dynamic air removal procedure<sup>2,3</sup> (with sufficient product drying)
- › steam sterilization according to EN 13060/EN 285 or ANSI AAMI S79 (for USA: FDA clearance)
- › validated according to EN ISO 17665 (valid IQ/OQ (commissioning) and product specific performance qualification (PQ))
- › maximum sterilization temperature 138 °C (280 °F; plus tolerance according to EN ISO 17665)
- › sterilization time (exposure time at the sterilization temperature):

Area	fractionated vacuum/dynamic air removal	gravity displacement
USA	at least 4 min at 132 °C (270 °F), drying time at least 20 min <sup>4</sup>	not recommended
other countries	at least 3 min <sup>5</sup> at 132 °C (270 °F) / 134 °C (273 °F), drying time at least 20 min <sup>4</sup>	not recommended

<sup>2</sup> at least three vacuum steps

<sup>3</sup> The less effective gravity displacement procedure must not be used in case of availability of the fractionated vacuum procedure, will require significantly longer sterilization times and is to be validated dependent on product, packaging, sterilizer, program, and parameters under sole responsibility of the user.

<sup>4</sup> The effectively required drying time depends directly on parameters in sole responsibility of the user (load configuration and density, sterilizer conditions, ...) and by this is to be determined by the user. Nevertheless, drying times less than 20 min must not be applied.

<sup>5</sup> respectively 18 min (inactivation of prions, not relevant for USA)

The fundamental suitability of the instruments for an effective steam sterilization was demonstrated by an independent, governmentally accredited and recognized (§ 15 (5) MPG) test laboratory by application of the steam sterilizer HST 6x6x6 (Zirbus technology GmbH, Bad Grund) and the fractionated vacuum/dynamic air removal procedure. For this, typical conditions in clinic and doctor's practice as well as the specified procedure were considered.

The flash/immediate use sterilization procedure must not be used. Do not use dry heat sterilization, radiation sterilization, formaldehyde and ethylene oxide sterilization, as well as plasma sterilization.

##### Storage

Please store the instruments after sterilization in the sterilization packagings at a dry and dust-free place.

##### Material resistance

Please take care that the listed substances are not ingredients of the cleaning or disinfection detergent:

- › organic, mineral, and oxidizing acids (minimum admitted pH-value 5.5)
- › strong lyes (maximum admitted pH-value 11, neutral/enzymatic or alkaline cleaner recommended)<sup>6</sup>
- › organic solvents (for example: acetone, ether, alcohol, benzine)
- › oxidizing agents (for example: hydrogen peroxide)
- › halogens (chlorine, iodine, bromine)
- › aromatic, halogenated hydrocarbons

<sup>6</sup> For the bur blocks alkaline cleaners must not be applied (maximum admitted pH-value 9).

Please do not clean any instruments and bur blocks by use of metal brushes or steel wool.

Please do not expose any instruments and bur blocks to temperatures higher than 142 °C (288 °F)!

Please do not apply acidic neutralizing agents or cleaning aids.

##### Reusability

The instruments can be re-used – in case of adequate care and if they are undamaged and clean as indicated in chapter "Specific aspects". The user is responsible for each further use as well as for the use of damaged and dirty instruments (no liability in case of disregard).

**Attachment A: Specific aspects**

- Diamond products and ceramic grinding tools:
  - Use particular care during the cleaning of the grinding surfaces and ensure that all residues are removed
- Bur blocks/instrument trays:
  - Cleaning and disinfecting only without products being loaded (products must not be cleaned and disinfected whilst they are in the bur block/ instrument tray)

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*The document in other languages as well as the latest and applicable version you find on the website [www.jota.ch](http://www.jota.ch).*

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Instrument group	brush	specific/additional procedure in case of				packing	sterilization	maximum admitted cycle number (confirmed by validation, but dependent on specific application)	recommended classification according to KRINKO/RKI/BfArM guidance (only German, with respect to intended use)
		pretreatment	manual cleaning/ disinfection	automated cleaning/ disinfection	maintenance				
stainless steel instruments	standard	standard	standard	standard	lubrication <u>not</u> admitted	standard	standard	10	critical B
regular steel instruments	standard	standard	standard	standard	lubrication <u>recommended</u>	standard	standard	10	critical B
silicone polisher	standard	standard	standard	standard	lubrication <u>not</u> admitted	standard	standard	5	critical B
endodontic instruments without stopper	endodontic brush	standard	standard	standard	lubrication <u>not</u> admitted	use of bur blocks not admitted	use of bur blocks not admitted	10	critical B
endodontic instruments with stopper	endodontic brush	mounted	mounted move the stopper at least three times during disinfection	mounted	lubrication <u>not</u> admitted	use of bur blocks not admitted	use of bur blocks not admitted	10	critical B
all other instruments	standard	standard	standard	standard	lubrication <u>not</u> admitted	standard	standard	10	critical B

**General application and safety instructions**

for the medical device from Jota AG  
Issued: January 2017

- Jota AG products (dental, maxillary surgery, general surgery,) must only be used by dentists, doctors and/or the respective experts who, due to their training and experience, are intensely familiar with the use of these products and who have the corresponding expertise in the respective specialist fields. The use of surgical products requires relevant expertise and experience in dental implantology, maxillary surgery and/or other surgical fields including diagnosis, preoperative planning and surgical techniques.
- It is the sole responsibility of the doctor in charge who, depending on the respective situation (indication), decides on the actual use of the Jota AG products for each individual case
- All JOTA AG products have been developed for specific applications. Therefore, inappropriate use can result in the premature wear and tear of the products and put patients and users at risk.

**Application**

- In order to avoid damaging the instruments, they must be removed from the blister pack by pulling off the back-sheet.
- It is essential to only use turbines as well as hand and angle pieces that are technically and hygienically faultless, maintained and cleaned.
- The instruments must be rotating when applied on material. They should not be placed on material and then brought to rotation.
- Rotating instruments need to be clamped as far down as possible with their speed set before applying them on the object. are used with the rotary instruments.
- Using the instruments for canting or levering should be avoided as it increases the risk of breakage.
- Depending on the application, it is recommended to use protective goggles while using the instruments. Users of diamond disks should use a disk protection device.
- Inappropriate use of the products leads to badly executed work and increased risk.
- When working with dry materials, it is recommended to use a suction cleaning device.
- In particular, users of hand tools should take care to use them gently and with consideration.
- The user must at all times avoid touching the instruments and parts without protection (protective gloves should be worn).
- Thermal bone damage caused by rotating and oscillating tools (e.g. pilot burr, conical burr, expansion burr) must at all times be avoided (user training, working at low speed and with sufficient cooling).
- During intraoral application attention has to be made to the fact that the products are protected against aspiration or falling on the floor.

**Use of pressure**

- Users of the instruments should at all times avoid applying excessive pressure. This can damage the working part of the instruments and cause the cutting edges to break off. At the same time, it generates excessive heat.

- The use of excessive pressure when using grinding tools can cause the abrasive particles to break off or the instrument to become clogged and lead to heat generation.
- During polishing, excess pressure can lead to heat generation.
- Due to overheating, excess pressure can damage the dental pulp or, due to broken off cutting edges, it can result in undesired rough surfaces. In such cases, even instrument breakage cannot be excluded.

**Cooling**

- In order to avoid excessive heat generation during preparation, a sterile water/sodium chloride solution supplied via a permanent external feeding device should be used to ensure sufficient cooling during use of the instruments.
- When using FG instruments that are more than 22 mm long or whose head diameter exceeds 2 mm, additional external cooling is required.
- Insufficient cooling will lead to irreversible damage to the bone and/or the adjacent tissue.

**Storage, disinfection, cleaning and sterilization**

- Unless there is explicit information to the contrary, all Jota AG products are supplied in non-sterile packaging and, depending on the application, they need to be sterilised prior to use. Prior to their first use on the patient and immediately after each use, all products need to be disinfected and sterilised. Inappropriate cleaning and sterilising of the instruments can result in the patient being infected with harmful bacteria.
- You will find detailed instructions for the disinfecting, cleaning and sterilising of products in the Instructions for the processing of instruments produced by Jota AG on the previous page. We would also be happy to provide you with these instructions at your request. They are also available on the internet at [www.jota.ch](http://www.jota.ch).
- The products should be stored in appropriate, hygienically maintained containers. The same applies to sterilised instruments. The stored products must be protected from dust, humidity and recontamination. Instructions as to maximum storage duration must be adhered to.

**Speed recommendations for rotary instruments**

- Following the instrument-specific speed recommendations produces the best results.
- Exceeding the maximum admissible speed (rpm) when using long and pointed instruments tends to produce vibrations that can lead to the destruction of the instrument.
- When using working parts with diameters exceeding the thickness of the shaft, excessive speed can release great centrifugal forces that may cause the shaft to bend and/or the instrument to break. Therefore, the maximum admissible rpm must never be exceeded.
- Please consult the manufacturer's information (see catalogue or [www.jota.ch](http://www.jota.ch)) for the recommended and the maximum admissible speed ranges. Non-compliance with the maximum admissible speed puts safety at risk.
- Generally, the following rules apply:

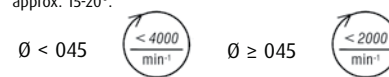
- The larger the working part of an instrument the lower the speed
- Surgical instruments: suitable for geared down micro-motor hand and angle pieces 10:1 with stable ball bearings. Speed 600 to 800 rpm with physical and, possibly, sterile external cooling or internal cooling when using the respective hand piece.

**Discarding worn instruments and parts**

- Jota AG products can principally be reused several times – unless specifically indicated and labelled otherwise. Rotating instruments are subject to wear. The option of and accountability for multiple use of a product and the frequency of application is solely the decision and own responsibility of the treating clinician based on the application in each case and the possible wear of the products. If in doubt, the products should always be sorted out early and replaced.
- Broken off cutting edges of instruments cause vibrations and great forces of pressure, which, in turn, leads to broken preparation corners and rough surfaces.
- Bare patches on diamond instruments indicate a lack of abrasive particles and can be a sign of blunt instruments. This leads to excessive temperatures during instrument use.
- Instruments that are bent and/or do not run true should be discarded forthwith.
- With the reuse of disposable products the risk of infection cannot be excluded and a risk-free functional safety cannot be guaranteed.

**Additional instructions regarding the use of trepans**

- When using trepans, you have to proceed with particular care. For example, it is advisable not to exceed the recommended rpm speed ranges.
- In order to prepare for the actual use of a trepan, it should be set to produce counter-clockwise rotations creating a groove in the bone. Afterwards the trepan can be inserted into this groove and, using clockwise rotations, it can be moved further down.
- Carrying out a prior X-Ray is essential to establish the maximum possible drilling depth and to maintain the necessary distance, for example, to the mandibular nerve. As an additional safety measure to spare the nerve, the axial direction of the trepan countersink attachment, based on the sagittal level of the ascending branch, must be milled laterally at an angle of approx. 15-20°.

**Further comments:**

- Due to statutory regulations, returned goods can, on principle, only be accepted if the complete batch number is provided. This number can be found on the product packaging.

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© JOTA AG, Switzerland All rights reserved. Reproduction, also by extract, are only permitted with written authorization of JOTA AG, Rüthi. All offers, orders and deliveries are subject to our „General Sales and Delivery Terms“. We reserve the right to modify our range of products and their design as well as to deviate from the illustrations and data shown.